

**National Institute for Occupational Safety and Health
Respirator Branch
Test Data Sheet**



Task Number: TN-25934

Reference No.: CFR 84.172

Test: Exhalation Resistance Test

STP No.: 3

Manufacturer: Sleepnet Corporation

Filter Type: Filter Only

Tested - #48012, #48013 Mask Body -Regular & X-Large
#39826 Envo Pro N95 Filter

Item Tested: 48012, 48013

Sample	Maximum Allowable Resistance (MM of H2O)	Actual Resistance (MM of H2O)	Result
	Exhalation	Exhalation	
1	25	9.7	PASS
2	25	9.1	PASS
3	25	9.7	PASS

Overall Result: PASS

Comments:

Sample 2 is the 48013 XL.

Was all equipment verified to be in calibration throughout all testing? Yes No

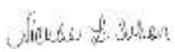
Signature:

Date: 8/9/2022

Engineering Technician

Task Number: TN-25934
Test: Exhalation Resistance Test
Manufacturer: Sleepnet Corporation

Reference No.: CFR 84.172
STP No.: 3

Signature: 

Date: 8/9/2022

Engineering Technician

**National Institute for Occupational Safety and Health
Respirator Branch
Test Data Sheet**



Task Number: TN-25934
Test: Inhalation Resistance Test
Manufacturer: Sleepnet Corporation
Item Tested: 48012, 48013

Reference No.: CFR 84.172
STP No.: 7

Filter Type: Filter Only

Sample	Maximum Allowable Resistance (MM of H2O)	Actual Resistance (MM of H2O)	Result
	Inhalation	Inhalation	
1	35	10.4	PASS
2	35	7.6	PASS
3	35	10.4	PASS

Overall Result: PASS

Signature:

Date: 8/9/2022

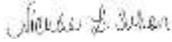
Engineering Technician

Task Number: TN-25934
Test: Inhalation Resistance Test
Manufacturer: Sleepnet Corporation
Item Tested: 48012, 48013

Reference No.: CFR 84.172
STP No.: 7

Comments:
Sample 2 is the 48013 XL.

Was all equipment verified to be in calibration throughout all testing?	<input checked="" type="radio"/> Yes	<input type="radio"/> No
---	--------------------------------------	--------------------------

Signature: 

Date: 8/9/2022

Engineering Technician